# Istruzioni d'uso

ver. 2 del 15/03/2021

Biologia Molecolare

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## CMV REAL TIME

<b>REF</b>	AA1489/48 GA028/2
	48 TEST 2 TEST
<b>CND</b>	W0105040205

**CE**  
0459

**IVD**

**NUCLEAR LASER MEDICINE S.r.l.**

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The screenshot displays the software's configuration steps. The first section, 'What type of experiment do you want to set up?', offers options like 'Quantitation - Standard Curve', 'Relative Standard Curve', 'Comparative Ct', 'Melt Curve', 'Genotyping', and 'Presence/Absence'. The second section, 'Which reagents do you want to use to detect the target sequence?', provides 'TaqMan Reagents', 'SYBR Green Reagents', and 'Other'. The third section, 'Which ramp speed do you want to use in the instrument run?', allows choosing between 'Standard (~ 2 hours to complete a run)' and 'Fast (~ 40 minutes to complete a run)'. Each section includes a brief explanatory note.

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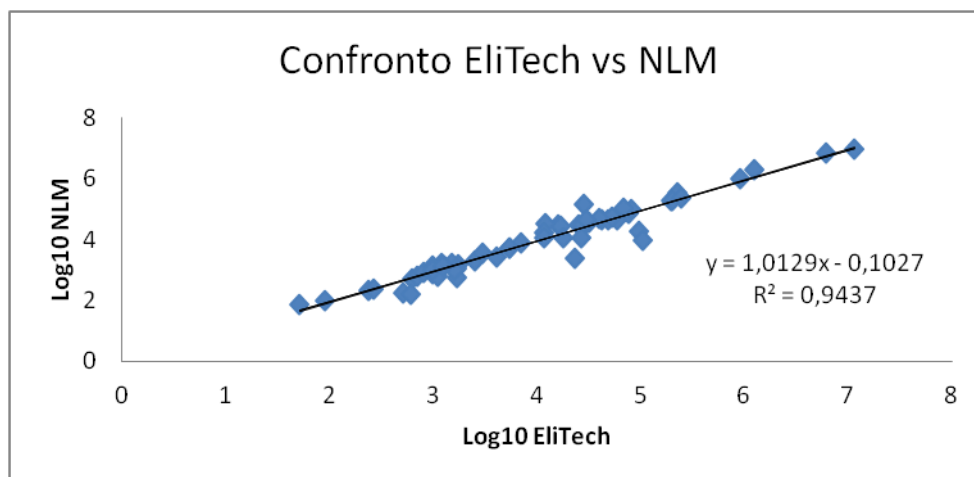
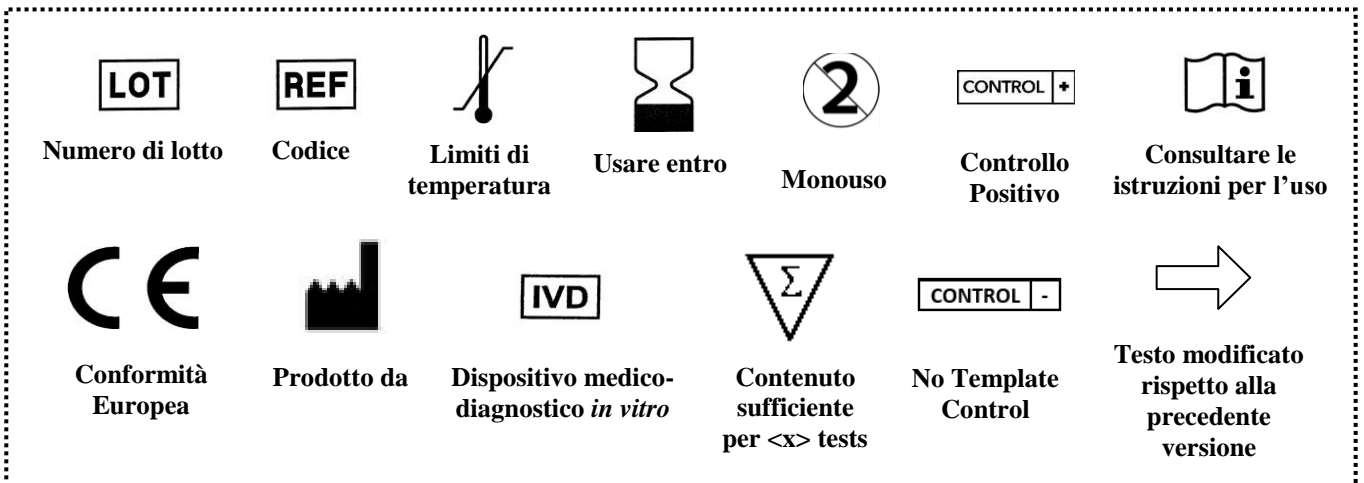
**IVD**

**NUCLEAR LASER MEDICINE S.r.l.**

HEADQUARTER: Via Cascina Conighetto – BUSINESS OFFICES: Viale delle Industrie, 3 - 20049 SETTALA (MI), ITALY

Phone: (+39) 02/95.24.51 - Fax (+39) 02/95.24.52.37 - (+39) 02/95.24.52.38 WEB: [www.nlm.it](http://www.nlm.it) - E-MAIL: [segreteria@nlm.it](mailto:segreteria@nlm.it)

ISO 9001 Quality Management System Organization certified and ISO 13485 Medical Sector Quality Management System certified (IMQ certification body - CSQ certification).

## ⇒ INTENDED USE

**CMV REAL TIME** provides reagents for the qualitative determination of the human Cytomegalovirus (CMV DNA) in amniotic fluid, urine and CSF specimens, and for qualitative and quantitative analysis of blood and plasma samples by Real Time PCR.

- Qualitative analysis involves the use of a Positive Control for the amplification step and of a Negative Control (NTC) that can be extracted together with the samples or directly amplified.
- Quantitative analysis involves the use of **CMV REAL TIME** kit together with **CMV Standard curve** kit (NLM code GA028/2), supplied separately.

This product is validated and compatible with the extraction kit on column (NLM code AA1001) and the most common automatic extraction systems: MagCore (NLM code AA1185 and AA1186), QiaSymphony (NLM codes AA1439/192 and AA1440/96).

The validated and compatible instruments for Real Time PCR are: CFX 96 (BioRad), ABI7500 Fast (Applied Biosystems) and RotorgeneQ (Qiagen).

**CMV REAL TIME** is indicated, together with the other laboratory parameters and the medical case of the patients, for the clinical management of CMV patients.

The test can be used both for a qualitative analysis and to evaluate the response to pharmacological treatment, expressed by the variation of the amount of viral DNA in biological samples.

The device includes the use of the enzyme Uracil-N-Glycosidase (UNG) to prevent contamination and false positive results.

The kit is to be considered for professional use only.

## INTRODUCTION

Cytomegalovirus (CMV) is a globally spread virus belonging to the Herpesvirus family. The virus is very common and can infect anyone. Once the infection is contracted, the virus remains latent within the body for life, but can reactivate if the immune system weakens. It can also cause mortality in immunocompromised patients, specifically in transplant recipients. The most important aspect related to CMV is congenital infections. Infection contracted during pregnancy and transmitted to the fetus can in fact cause serious damage to the child. Man is the only reservoir of CMV infection, whose transmission occurs from person to person through body fluids, including blood, saliva, urine, seminal fluids, vaginal secretions and milk. The infection can occur by contact person-person, by mother-fetus transmission during pregnancy or mother-child during lactation, for transfusions and transplants of infected organs. HCMV has a prolonged replication cycle; it has the most extensive genetic heritage (236 kbp) of all known human viruses. The genome is linear with a double-stranded DNA molecule. Viral load is routinely analyzed for diagnosis and choice of therapy.

## REFERENCES

1. Cunningham et al., **Sequences of complete human cytomegalovirus genomes from infected cell cultures and clinical specimens**, Journal of General Virology (2010), 91, 605–615
2. Jong Eun Park et al., **Performance Evaluation of the Real-Q Cytomegalovirus (CMV) Quantification Kit Using Two Real-Time PCR Systems for Quantifying CMV DNA in Whole Blood**, Ann Lab Med 2016;36:603-606
3. Sara Jones, Erika M. Webb, Catherine P. Barry, Won S. Choi, Klara B. Abravaya, George J. Schneider, Shiaolan Y. Ho, **Commutability of Cytomegalovirus WHO International Standard in Different Matrices**, Journal of Clinical Microbiology June 2016 Volume 54 Number 6
4. Camille N. Kotton et al., **Updated International Consensus Guidelines on the Management of Cytomegalovirus in Solid-Organ Transplantation**, Transplantation Volume 96, Number 4, August 27, 2013
5. Pillet et al., **Comparative evaluation of the QIASymphony RGQ system with the easyMAG/R-gene combination for the quantitation of cytomegalovirus DNA load in whole blood**, Virology Journal 2012, 9:231

## ⇒ PRINCIPLES OF THE PROCEDURE

The test is based on Real-Time PCR technique that includes two steps:

1. Viral DNA Extraction
2. Amplification and detection of the target sequence by Real Time PCR

The Internal Control (GAPDH DNA) is purified with the target and allows to monitor the assay. In case of poorly cellularised matrix (urine, amniotic fluid, CSF) it is possible to add the exogenous Internal Control directly to the sample before the extraction step.

Two different controls are supplied: a negative control, for the extraction and amplification steps, and a positive control specific only for the PCR phase.

For quantitative analysis a panel of four calibrators is provided in liquid ready-to-use form, supplied separately (NLM code GA028/2).

### 1. DNA Extraction

⇒ The DNA of CMV and Internal Control can be extracted by different kind of specimens (blood, plasma, urine, cerebrospinal fluid and amniotic fluid), using manual extraction through column system (NLM code AA1001) or by automatic system like MagCore (NLM code AA1185 and AA1186,) and QIASymphony (NLM code AA1439/192 and AA1440/96).

### ⇒ 2. Amplification and Detection

#### Specific target for the amplification

The amplification and the detection of the **MIEA gene of CMV** (major immediate early antigen, HCMVUL123) allow to determine the viral DNA presence (amplicon: 133 bp).

At the same time the device amplifies the intron 2 of GAPDH gene (amplicon: 234 bp) as endogenous or exogenous Internal Control.

The Internal control is extracted and amplified in each sample to monitor the whole procedure.


The positive control and calibrators consist of cloned synthetic DNA of MIEA UL123 gene and are supplied in liquid ready-to-use form.

#### Amplification


Following the extraction step DNA of CMV and IC is amplified by Real-Time PCR, using appropriate enzymes and specific reagents.

#### Detection

The test is based on Real Time PCR technique that allows to monitor in real time the amplification of the target sequences using appropriate probes marked with fluorescent molecules. In the amplification mix there are two specific probes respectively for the viral genome target region and for a human genome region that acts as Internal Control (GAPDH), each labeled with a donor fluorophore and a fluorophore acceptor (according to the "dual labeled" chemistry probe). The donor fluorophore is different for the two probes: FAM for CMV and JOE for IC. In the presence of a light source, when the probes are intact the fluorescence emitted by the donor fluorophore is absorbed by the acceptor fluorophore. During amplification, the probes pair up each to their specific target and subsequently are degraded by the 5' → 3' exonuclease activity of the DNA polymerase. When the donor and acceptor fluorophores are separated, the donor fluorescence can be detected at its specific wavelength: in this way the amplification of the viral DNA and IC can be monitored during the reaction.

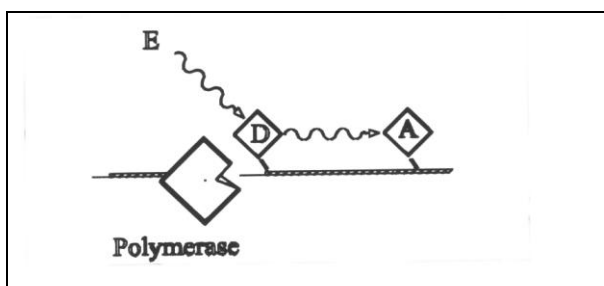a) Energy (E) emitted from the "donor" fluorophore (D) is absorbed by the "acceptor" (A) fluorophore.

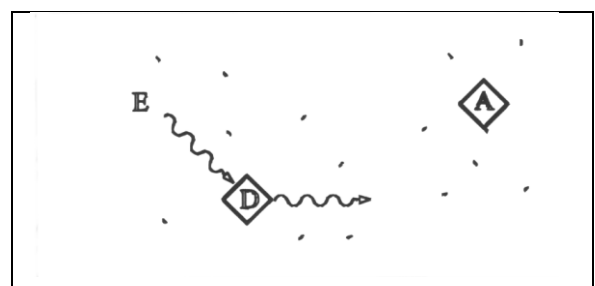b) The exonuclease activity of the DNA polymerase separates the fluorophore D from the fluorophore A through hydrolysis with a consequent increase in the fluorescence signal.

## ⇒ REAGENTS COMPOSITION (Store at -25/-15°C)

<b>AA1489/48</b>	<b>Quantity</b>
<b>CMV Mix</b> (including primers, probes, amplification enzyme and UNG)	2 x 500 µl
<b>Positive Control (CP CMV)</b> (synthetic DNA of target region; liquid ready-to-use form)	2 x 30 µl
<b>No template Control (NTC)</b> (Water to be used as NTC from extraction step)	2 x 700 µl
<b>Exogenous Internal Control (exIC)</b> (Synthetic DNA of IC; liquid ready-to-use form)	2 x 350 µl

<b>GA028/2</b>	<b>Quantity</b>
<b>CMV Standard curve</b>	<b>2 set</b>
CALIBRATORE 1 CMV (Synthetic DNA, liquid ready-to-use form)	2 x 30 µl
CALIBRATORE 2 CMV (Synthetic DNA, liquid ready-to-use form)	2 x 30 µl
CALIBRATORE 3 CMV (Synthetic DNA, liquid ready-to-use form)	2 x 30 µl
CALIBRATORE 4 CMV (Synthetic DNA, liquid ready-to-use form)	2 x 30 µl

The components volumes show above refer to the standard size of the kit.

Reduced packaging of the kit are available on request for evaluation and/demonstration of the product.

## ⇒ STABILITY AND STORAGE

- CMV Mix, NTC, Exogenous Internal Control, CP and Calibrators are stable until the expiration date on the label when stored at the indicated temperature (-25/-15° C).
- Thaw all reagents on ice or at +2/+8° C.
- Mix can be subjected to maximum of five freeze-thawing cycles. Prepare aliquots if more freeze-thawing cycles are expected.
- The probes in the CMV Mix are photosensitive: avoid prolonged exposure to light.
- Positive control and calibrators are ready to use; it is recommended to dispense them into the tubes / amplification plate after the samples, to avoid possible contamination. Each aliquot is disposable: after the first use discard the remaining volume.

## ⇒ PRECAUTIONS

- Only professional and opportunely trained personnel should carry out the procedure using good laboratory practices and common protective equipment.
  - The use of the following solutions, together with the presence of UNG enzyme, allows to minimize the risk of cross-contamination:
    - ✓ It is recommended to perform the assay in two different areas:
      - Area 1: pre-PCR (samples handling, extraction and PCR set up)  
Pay attention not to contaminate the reagents during handling.
      - Area 2: post-PCR (Real Time PCR)
- Each area should be provided with dedicated equipment and consumables (lab coats, centrifuge, tubes, pipettes, etc).
- ✓ Workspace should be organized to ensure that the workflow occurs in one direction, from pre-PCR area to post-PCR area.
  - ✓ At the end of the procedure it's recommended to clean the working surfaces (bench and vertical downflow airbox) with 5-10% bleach (final concentration of sodium hypochlorite: 0,5% w/v), freshly prepared. Follow the manufacturer recommendations to clean the instruments.

- ✓ All disposable items (tips and tubes) must be DNase, RNase free. Use aerosol-resistant pipette tips to avoid pipettes contamination. Use a new tip every time a volume is dispensed.
  - ✓ Use vertical downflow airbox with UV lamp
  - ✓ Change gloves frequently
  - ✓ The absence of contamination can be monitored by the use of the NTC in each run, starting from the purification step.
- Discard all used material in accordance with local and national regulations in force.
  - Do not eat, drink, smoke or apply cosmetics in areas where reagents or specimens are handled.
  - In case of contact with reagents rinse immediately with water and seek medical advice.
  - Do not use device after its expiration date.
  - Do not mix reagents from different lots.
  - Do not use the device if the box is damaged; contact the supplier.
  - It is advisable to have constant and uniform laboratory temperature, avoid to place the instruments near heating/cooling sources that may compromise the right working.

## ⇒ MATERIAL REQUIRED BUT NOT PROVIDED

### AREA 1

DNA extraction kit  
 Vertical downflow airbox  
 Dedicated adjustable volume pipettes set and aerosol barrier tips  
 Filter tips  
 2 ml and 1,5 ml DNase, RNase free tubes

### AREA 2

CMV Standard Curve (NLM code GA028/2) for quantitative analysis  
 Vertical downflow airbox  
 Dedicated adjustable volume pipettes set  
 Aerosol barrier tips DNase, RNase free  
 DNase, RNase free tubes, strips or plates, specific for the Real Time PCR instrument used  
 Real Time PCR Instrument

## COLLECTION AND STORAGE

- Whole blood and plasma samples must be prepared and stored according to laboratory directions.  
 Unless otherwise stated, it is suggested to store at +2/+8° C for up to three days and do not freeze them to avoid cell lysis and loss of bacterial DNA titre. It is also advised not to freeze for plasma separation from whole blood samples.
- Urine samples: samples of urine from the first jet must be collected, in containers without preservatives as indicated by the laboratory, transported at room temperature (+15/+25° C) for up to four hours or stored at +2/+8° C for a maximum of three days. If possible, avoid freezing of urine samples of the first jet. Freezing can cause the precipitation of inhibitors, cell lysis and, consequently, the loss of bacterial DNA
- Cerebrospinal and amniotic fluid samples: they must be prepared and stored according to the the laboratory instruction. Store at -25/-15° C.

## ⇒ ASSAY PROCEDURE

### VIRAL DNA EXTRACTION

#### Extraction methods:

- MANUAL EXTRACTION **cod. NLM AA1001.**
- AUTOMATIC EXTRACTION **cod. NLM AA1185, AA1186, AA1439/192 and AA1440/96.**

For preparation, use and disposal of the reagents, refer to the specific instructions for use. The possible specimen to use are: EDTA whole blood, plasma, amniotic fluid, cerebrospinal fluid, urine. For samples preparation it is necessary to follow the protocol concerning the specific method.

For beginning and elution volumes refer to table below:

Extraction method	Starting volume	Elution volume
Manual extraction (code NLM <b>AA1001</b> ).	200 µl	100 µl
Automatic extraction with MagCore (Diatech/MGB 400-02, code NLM <b>AA1185</b> and Diatech/ MVN400-04, code NLM <b>AA1186</b> )	400 µl	100 µl
Automatic extraction with QIAasympyphony (Qiagen/937055, code NLM <b>AA1440/96</b> – protocol <b>Cell free 500</b> )	700 µl	110 µl
Automatic extraction with QIAasympyphony (Qiagen/937236,-code NLM <b>AA1439/192</b> – protocol <b>VirusBlood200 DSP</b> )	300 µl	

Store the extracted DNA at -25/-15° C, if it is not used immediately for amplification.

In case of extraction from poorly cellularised matrix (urine, amniotic liquid, CSF) it is possible to add the Exogenpus Internal Control (ex IC), supplied with the device.

Add 15 µl of ex IC to biological matrix with every extraction method independently from the sample volume in order to ensure the internal control amplification. In case of CMV positive sample the addition of ex IC does not affect the specific target signal.

## ⇒ REAL-TIME PCR

- Set the thermal profile before preparing the mix.
- **Mix gently the CMV Mix before using it.**
- Dispense 15 µl of **CMV Mix** for the number of samples to analyze and for Positive Control (for qualitative analysis) or the four calibrators (for quantitative analysis) and for NTC.
- Add 15 µl of DNA, of NTC, of PC or calibrators to the CMV Mix, and mix carefully by pipetting.

In case of QIAasympyphony PCR setup, centrifuge the tubes before placing into the real time instrument to collect all drops. The solution has to be placed onto the bottom of the tube.

**CAUTION:** mix and dispense carefully samples and CMV Mix, avoiding the formation of bubbles. If possible, briefly centrifuge the plate or strips before placing them into the Real Time instrument. It is suggested to dispense calibrators and the positive control after the samples, to avoid cross contamination.


## RUN SET UP AND DATA ANALYSIS

### ROTORGENE Q

- ⇒ Click the button "Edit Samples" to set the position of the tubes and the "Type": Standard, Positive Control, NTC or Unknown.
- In the "Setting" function, set the reaction volume (30 µl) and the rotor used (36/72 wells); check that the Green and Yellow channels are shown in the "Channels" window, otherwise set them by clicking on "Create new".
- In the "View → gain optimisation" function set the autocalibration to 58°C for the above channels (**Min reading 5 FI; Max reading 10 FI**) and select "Perform optimisation before first acquisition".
- Set the following thermal profile in the "Profile" function:

Cycle	Cycle Point
45°C, 5 min	
95°C, 3 min	
Cycling (5 repeats)	95°C, 5 sec
	58°C, 20 sec
	72°C, 15 sec
Cycling (45 repeats)	95°C, 5 sec
	<b>58°C, 30 sec; acquiring to Green and Yellow</b>
	72°C, 15 sec

- Place the tubes in the rotor of the RotorGene and close the lid of the instrument.
- Start the experiment by clicking on "Start".
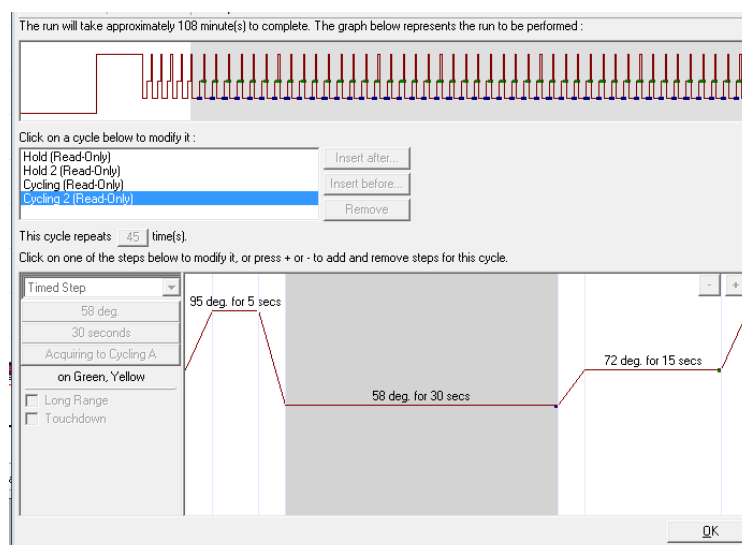
## DATA ANALYSIS

### Green Channel (CMV)


- Select "Analysis", "Quantitation", "Cycling A.Green", "Show".
- When the analysis window appears with the graph, select "Dynamic Tube" and "Slope Correct" and, if necessary, set "outlier removal" at 10%.
- Enter the value 0.1 in the "Threshold" window and, if necessary, delete the first 10-12 cycles: the intersection between the threshold line and the sample curve represents the Ct value.

Qualitative Analysis: the "Results" window shows the Ct values of the CMV positive samples and of the Positive Control

Quantitative Analysis: to set Standard Curve points select "Type" → "Samples" → "Standard". Then set the concentration for each calibrator, select "Given Conc." → "Samples" and enter the number of IU/ml (to choose the unit of measure go to "Unit" at the top right and choose "IU/ml" from the drop-down menu).

The screen of the standard curve will appear and the viral load will be assigned to each positive sample.


Example:
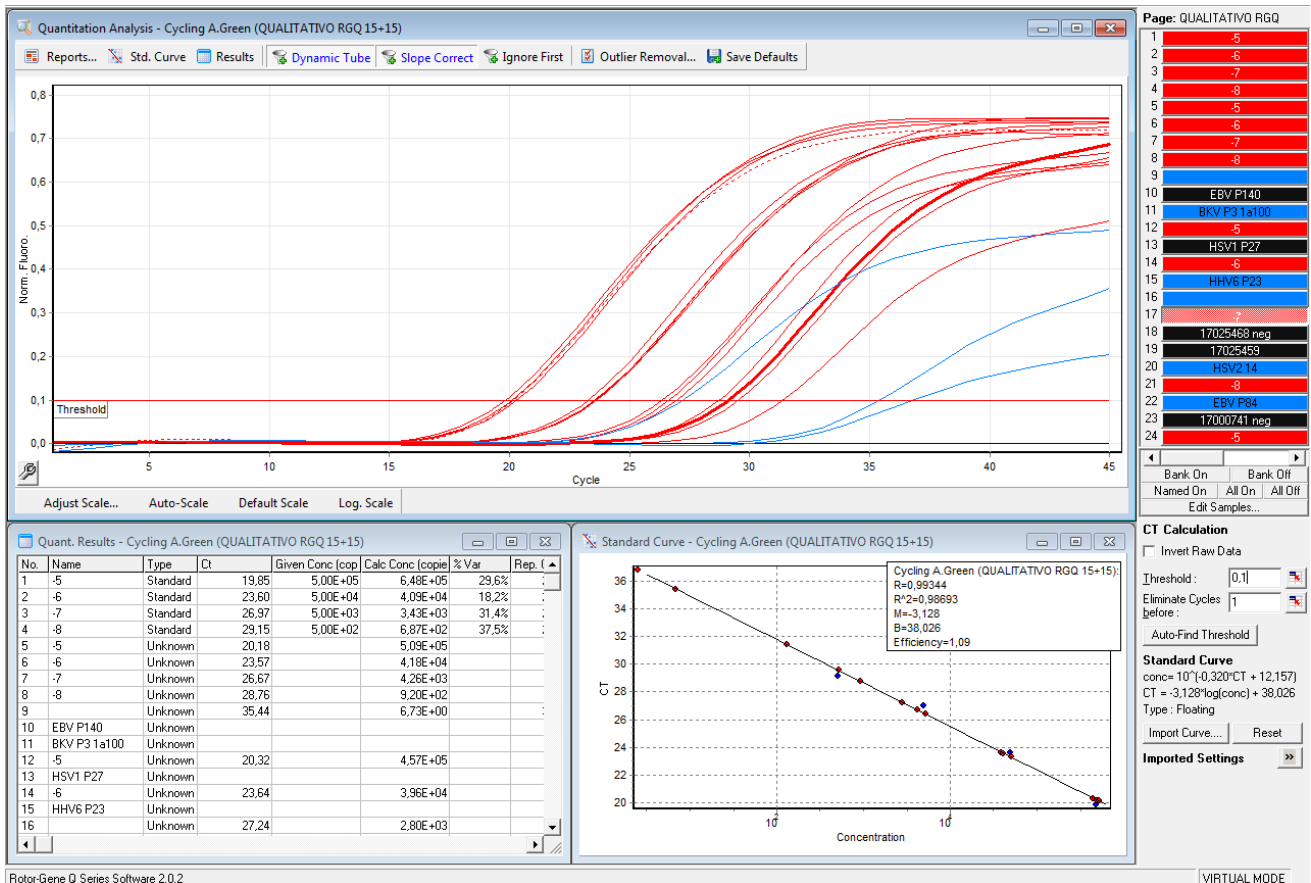
## Yellow Channel (IC)

- Select "Analysis", "Quantitation", "Cycling A.Yellow", "Show".
- When the analysis window appears with the graph, select "Dynamic Tube" and "Slope Correct", if necessary, set "outlier removal" at 10%.
- Enter the value 0.035 in the "Threshold" window and, if necessary, delete the first 10-12 cycles: the intersection between the threshold line and the sample curve represents the CT value.

Example:
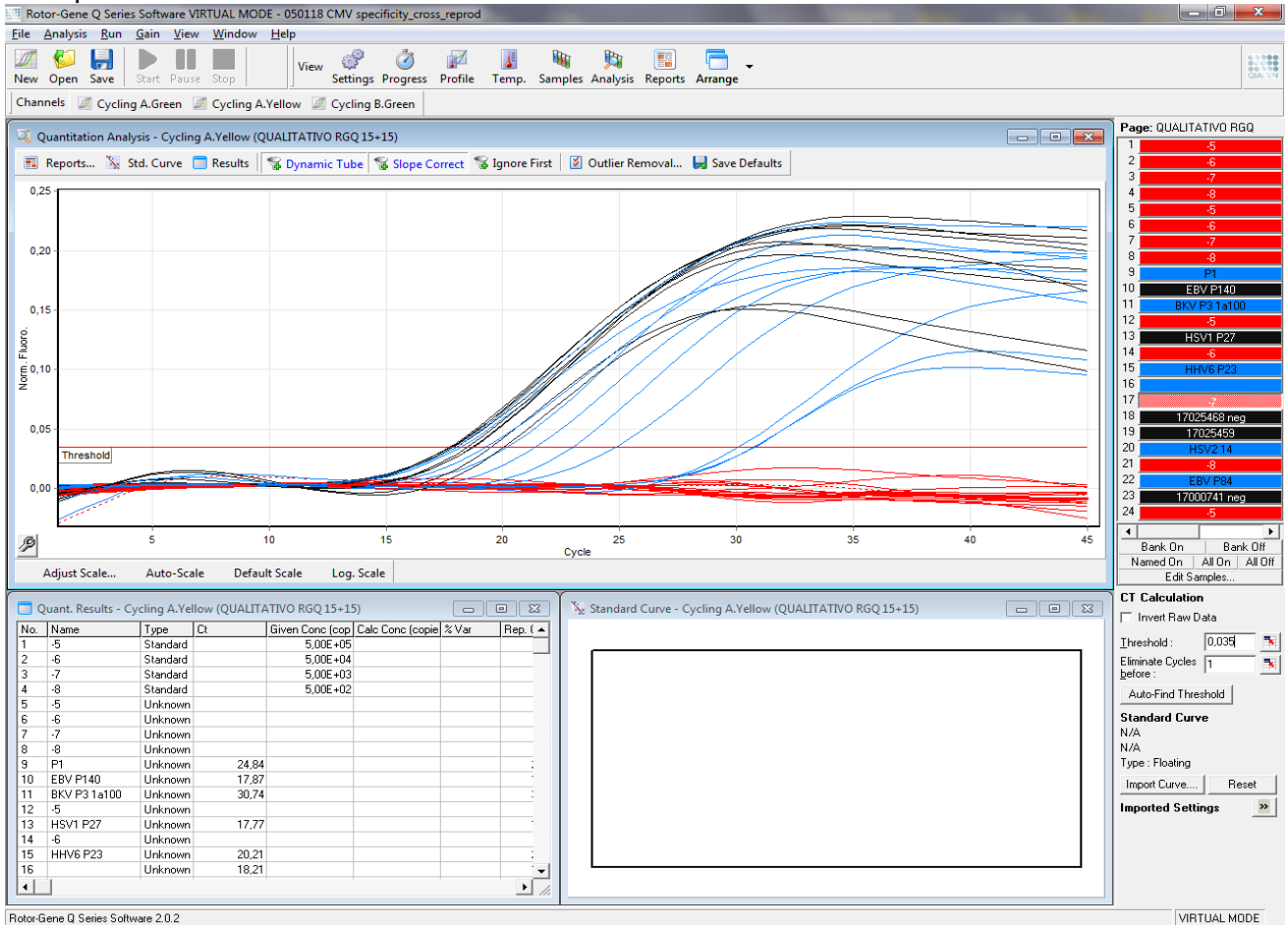
## CFX 96

- Open the Bio-Rad CFX Manager specific software.
- Select *Create new run* and choose the CFX model.
- Click OK.
- The Run Setup page opens

## PROTOCOL

It is possible to create a new protocol or select/modify an existing protocol.

### Create new protocol


- *Protocol* → *create New*
- Open the window *Protocol Editor – new*
- Insert the volume of the reaction in the *Sample Volume* box (30  $\mu$ l)
- Set the thermal protocol as below:

Cycle	Cycle Point
45°C, 5 min	
95°C, 3 min	
Cycling (5 repeats)	95°C, 5 sec
	58°C, 20 sec
	72°C, 15 sec
Cycling (45 repeats)	95°C, 5 sec
	<b>58°C, 30 sec; Plate Read</b>
	72°C, 15 sec

- To add a step, click *Insert Step* on the left side of the window. The step will be added to the position following the selected step.
- It is possible to change temperatures and time of each step, directly, from the screen or in the text below, by a double click on them.
- Check the presence of the symbol of a camera in the annealing step of the screen and the presence in the text below of the written “+ Plate Read”. If not, add them by clicking “Add Plate Read to Step”.
- Click OK. Save the new protocol.

### Select / modify an existing protocol

- It is possible to import a thermal profile from saved profiles stored in the PC (\* .pcrd): click *Select Existing* in the *Protocol page, Run Setup*.
- If you want to modify an existing protocol, click *Select Existing and Edit Select*. Then proceed as described above.

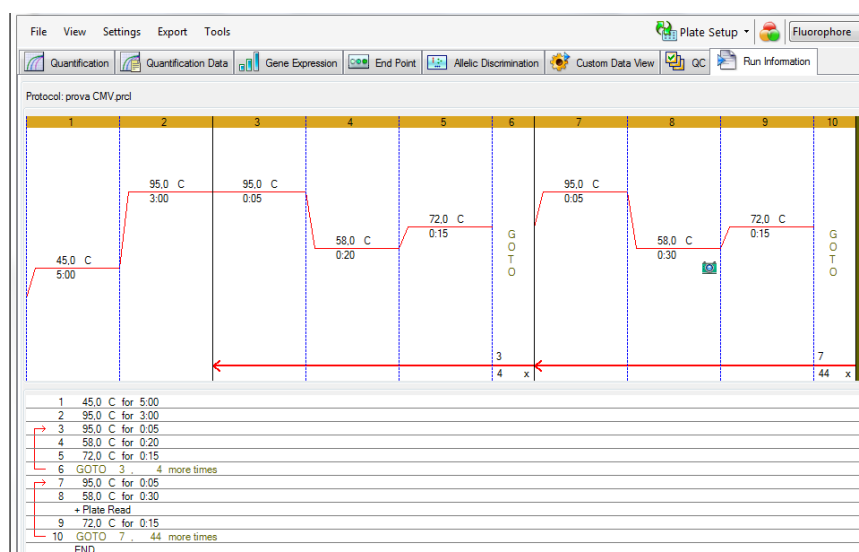
## PLATE

Select *Plate* in the window *Run Setup*.

It is possible to create a new plate or select/modify an existing plate.

### Create new plate

- *Plate* → *create New*
- Open the window *Plate Editor – new*
- Set the type of plate selecting *Settings* → *View/Edit Plate* → *Settings* → *Plate Type* and check *BR Clear*
- Select *Scan Mode* → *All channels*
- Click *Select Fluorophores* and choose FAM and JOE channels. (CAUTION: If JOE is not present, select the equivalent HEX fluorophore)
- Select the wells
- ⇒ • Set *Sample Type* → *Unknown, Standard, Positive Control or NTC* depending on the type of sample
- Write in *Target Name* the fluorophore name and click *load* to link the fluorophores to the wells: CMV for FAM and GAPDH for JOE/HEX. All samples will be analyzed for both fluorophores while Positive Control and calibrators will be analyzed only for FAM.
- Write the name of the samples in *Sample Name* and click *load* to associate the name with the wells.
- For the calibrators, enter the declared concentration values found in the box label and click on *Load* to confirm.
- Click *OK*. Save the plate.
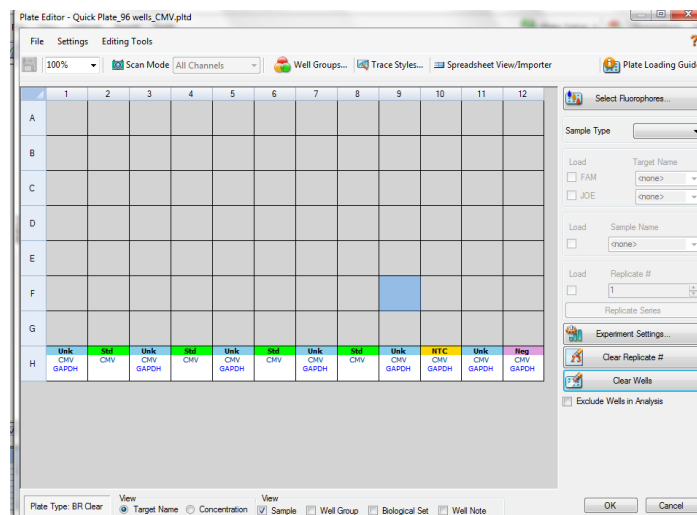
### Select/modify an existing plate

- It is possible to import a plate from the profiles saved and stored in the PC (\*.pcrd): click *Select Existing* on the page *Plate*.
- Click *Edit Select*. Then proceed as described above.

## RUN

- Select *Start Run* in the *Run Setup* window.
- The instrument is now ready to work. Prepare the plates / strips with CMV Mix and samples.
- Open the lid by clicking the *Open Lid* button and place the plate / strips in the instrument; close the lid by clicking the *Close Lid* button.
- Click *Start Run* to start the session.
- The *Run Details* window opens. Move to the *Real Time Status* page, to monitor the run.


## DATA ANALYSIS

At the end of the run the window *Data Analysis*→*Quantification* is shown.


To exclude a sample/calibrator from the analysis, go to “*plate setup-view/edit plate*” (at the top right of the screen), select the corresponding well and click on “*Exclude well in analysis*” at the bottom right.

At the top left of the screen select *Settings* and set:

- *Cq determination mode: single threshold*
- *Settings-Analysis Mode (version 1.6)/Baseline Setting (versions 3.x): Baseline subtracted curve fit* (for CFX Manager v1.6 and following versions)

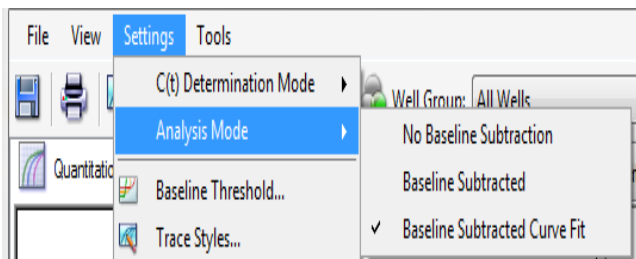
a) CFX Software v1.6

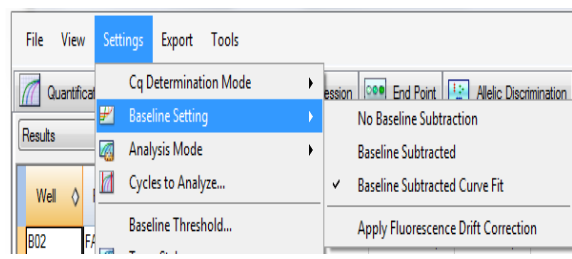
b) CFX Software v3.x

To set the threshold select one fluorophore at a time.

- **FAM:** *Baseline threshold: Single Threshold-User defined-set 200*. Click *OK*
- **JOE/HEX:** *Baseline threshold: Single Threshold-User defined-set 100*. Click *OK*
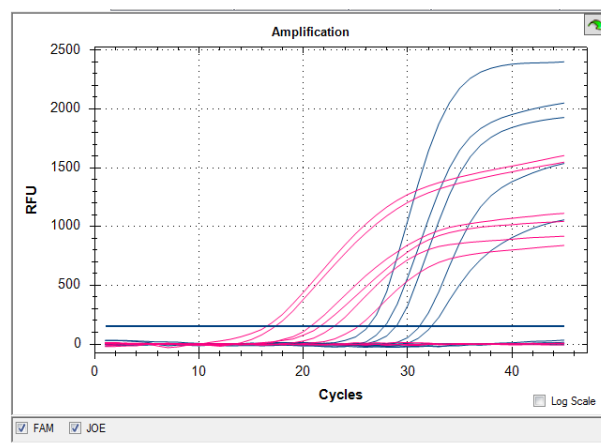
## Quantitative Analysis

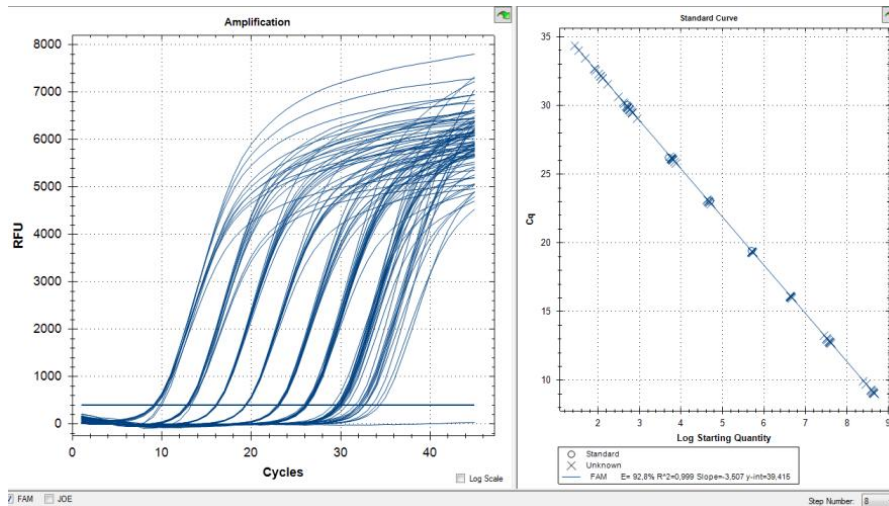
The screen on the left shows the amplification data (linear scale or log), the screen on the right shows the standard curve.

- To create a final report select “*Tools → Reports*” and select/deselect the information of interest.

**WARNINGS:** For a correct interpretation of the results, always pay attention to the screen of the run.

Example: in the image below the sample highlighted in green seems to show a positive Ct, but from the graphical data an instrumental artifact is evident. The result cannot be considered as valid and the sample must be repeated.
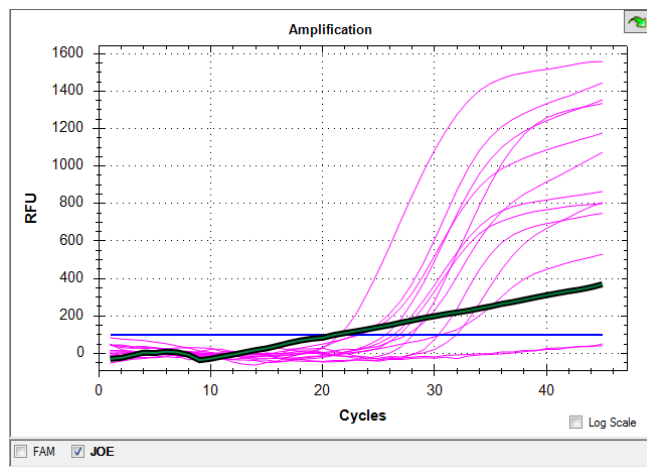
## ABI 7500 FAST

- Open the software (version v2.x) and click *Advanced Setup* on the left side of the screen to open a new window.
- In the Setup menu, on the left, select *Experiment Properties*, enter the experiment name and select the *7500 Fast (96 wells)*, *Quantitation Standard Curve*, *TaqMan® Reagents* and *Standard (2 hours to complete a run)*.

The screenshot displays three configuration panels in the software interface:

- What type of experiment do you want to set up?**
  - Buttons: Quantitation - Standard Curve (selected), Quantitation - Relative Standard Curve, Quantitation - Comparative Ct ( $\Delta\Delta Ct$ ), Melt Curve, Genotyping, Presence/Absence.
  - Text: Use standards to determine the absolute quantity of target nucleic acid sequence in samples.
- Which reagents do you want to use to detect the target sequence?**
  - Buttons: TaqMan® Reagents (selected), SYBR® Green Reagents, Other.
  - Text: The PCR reactions contain primers designed to amplify the target sequence and a TaqMan® probe designed to detect amplification of the target sequence.
- Which ramp speed do you want to use in the instrument run?**
  - Buttons: Standard (~ 2 hours to complete a run) (selected), Fast (~ 40 minutes to complete a run).
  - Text: For optimal results with the standard ramp speed, Applied Biosystems recommends using standard reagents for your PCR reactions.

## PLATE SETTINGS


- In the *Setup* menu, on the left, select *Plate Setup*
- Select *define targets and samples* at the top.
- In the left part of the page (define targets) select *add new target* and enter the name of the target; from the dropdown menu select FAM (for CMV) as a reporter and NFQ-MGB (non-fluorescence quencher-minor groove binder) as a quencher and a *color*. Repeat the operation inserting JOE/HEX as a reporter (for GAPDH).
- On the right side of the window (define samples) select *add new sample* for each sample to be analyzed and enter the *sample name*.
- On the top select *assign targets and samples*.
- Select the positions used in the plate and associate the fluorophore, checking the target on the left (assign). Also indicate in the *Task* entry the type of sample present in that particular position: *Positive Control*, *NTC* or *Unknown* (for samples), or, for quantitative analysis, consider the position for the *Standards*.
- For quantitative analysis: select the wells for the calibrators and click "S" in the "Task" window, then set the concentration in the "Quantity" box on the right. By clicking on Analyze, the software will assign a value expressed in IU/ml.
- Select *NONE* relatively *the dye to use as the passive reference*.

## Thermal Profile

- Select on the left *Run Method*.
- Select on the top *Graphical View*.
- Set 30  $\mu$ l in *Reaction volume for well*.
- Set the thermal profile
- Select *Add stage* Holding or Cycling from the drop-down menu to add a stage, while to add Step in the same stage use the Add Step menu.
- Ramps: All ramps must be left at 100%.
- Set fluorescence signal acquisition during annealing only in PCR cycles.
- Leave the *Enable AutoDelta* unselected

Cycle	Cycle Point
45°C, 5 min	
95°C, 3 min	
Cycling (5 repeats)	95°C, 5 sec
	58°C, 20 sec
	72°C, 15 sec
Cycling (45 repeats)	95°C, 5 sec
	<b>58°C, 30 sec; Plate Read*</b>
	72°C, 15 sec

\*Data collection on  (activated icon)

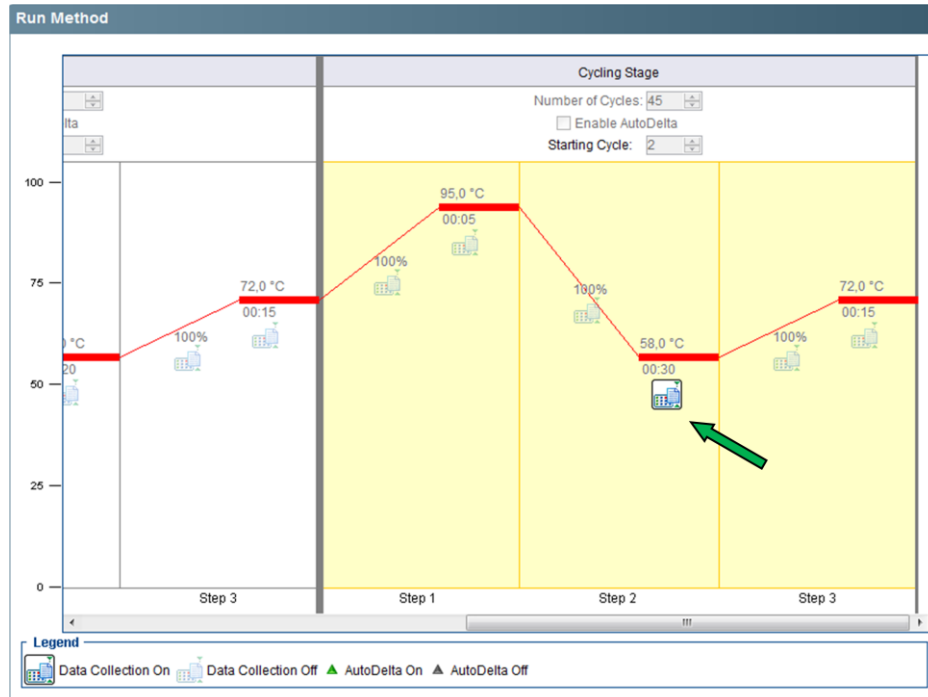- Check in the Run window that the Instrument Status is connected: the instrument is then ready to work.

## START RUN

- Prepare the plate with the CMV Mix and samples, seal it carefully taking care not to create air bubbles. Centrifuge the plate a few seconds to remove any bubbles.
- Open the instrument and place the plate on the support, with the A1 position at the top left.
- For any other indication, refer to the instrument user manual
- Save the run.
- To save the template, select File, Save as template. In the following sessions it will be possible to recall and modify the template by selecting File, New experiment, From template
- Press Start.

## DATA ANALYSIS

- Select Analysis on the left of the window
- To visualize amplification curves, select *Amplification plot* and set *dRn vs Cycle* for plot type, *Linear* for graph type and the *target* that you want to analyze.
- If you want to view one target at a time, choose the target of interest in “Options”.
- To modify Baseline and Threshold, select on the top right *Analysis, Analysis Settings* and open the Ct Settings window: select the target you want to modify:
  - o de-select *Use Default Settings*
  - o set these parameters:

	Plate	Strip
<b>FAM</b>	200.000	100.000
<b>JOE</b>	30.000	15.000
<b>Autobaseline selection</b>	<b>NO (Start=3; End=15)</b>	<b>NO (Start=3; End=15)</b>

- o select Apply Analysis Settings.
- click *Re-analyze*
- In the window *View Well Table* it is possible to visualize the Ct value and the viral load of each sample.
- Clicking on *Standard Curve*, on the left, it is possible to visualize the reaction efficiency, in case of Quantitative Analysis; otherwise it is possible to evaluate the run on the basis of the Ct of the samples and of the Positive Control (qualitative analysis).
- To generate a report, click on “Print report” and select the information of interest.

## ⇒ INTERPRETATION OF THE RESULTS

If the No Template Control (NTC) is analyzed in the run, interpret the results as follows:

	FAM Ct	JOE/HEX Ct	Result
<b>No Template Control</b>	Absent	Absent	Valid
	Absent	Present	Not Valid*
	Present	Present	Not Valid *

(\*)Hypothesis of contamination

## QUALITATIVE ANALYSIS

To establish the correct amplification it is necessary to evaluate the Positive Control signal, interpret the results as follows:

	Fam Ct	Joe/Hex Ct	Interpretation
<b>Positive Control</b>	Present	Absent	Valid result
	Absent	Absent	Not valid result*

(\*) The absence of FAM signal in the Positive Control could be due to DNA target degradation and/or procedure error.

For each type of samples perform the analysis as follows:

	Fam Ct	Joe/Hex Ct	Interpretation
<b>Sample</b>	Present	Present	CMV detected- valid result
	Absent	Present	CMV not detected – valid result
	Present/Absent	Absent	Repeat the test**

(\*\*) In case of poorly cellularized sample, add 15 µl of ex IC before the extraction phase.

## QUANTITATIVE ANALYSIS

- To proceed with the quantitative analysis it is necessary to evaluate the calibrators signals as follows:

	<b>Fam Ct</b>	<b>Joe/Hex Ct</b>	<b>Interpretation</b>
<b>Calibrators</b>	Present	Absent	Valid result
	Absent	Absent	Not valid result*

(\*) The absence of the FAM signal in the calibrators could be due to DNA target degradation and/or procedure error.

- To perform the quantitative analysis, verify that the Ct of the calibrators are included within the range showed on the box label (GA028/2) and assign the relative quantification values.
- If the interpolation line does not fit well with one of the four calibrators (eg in the presence of a Ct outside the corresponding acceptable RANGE shown inside the package), that calibrator can be excluded from the analysis, in a way to obtain a better approximation of the straight line to the remaining points; this allows a more precise quantification of the samples, without altering the clinical significance of the result.
- If the interpolation line is not suitable for more than one calibrator, the standard curve cannot be used and the run must be repeated.
- The concentration values out of the linear range of the assay cannot be considered accurate, so they must be interpreted with caution. If the calculated concentration is above the upper detection limit, report the result as CMV DNA > 5E + 08 IU / ml. To obtain a more accurate result, the test must be repeated using the starting sample diluted with a negative sample of the same matrix. Multiplying by the dilution factor applied, the initial concentration will be obtained.
- In case of samples with viral load lower than the limit of detection, the reproducibility of the results is not guaranteed.

## **TROUBLESHOOTING**

- In case of contamination it is important to recognize the source (the extraction and/or PCR phase). The use of the controls can be a valid aid to ensure the correct performance of the test and to identify the origin of this problem.
- Possible mutations in the target regions of the viral DNA can alter the primers and/or probes pairing, preventing the right identification of the sample.

## PERFORMANCE CHARACTERISTICS

### Diagnostic specificity

The diagnostic specificity of the kit CMV REAL TIME was determined by analyzing, following the instruction for use, 197 CMV negative samples of blood, plasma, urine, cerebrospinal and amniotic fluids. In agreement with these data, the diagnostic specificity is 99,5%.

### Analytical sensitivity

#### 1. Whole Blood

The analytical sensitivity was determined analyzing five scalar dilutions of a positive CMV blood sample, calibrated with the First International CMV WHO Standard (NIBSC code 09/162). For serial dilutions a negative CMV human blood was used. The analysis was carried out according to user manual, extraction included. The detection limits for CMV REAL TIME with whole blood samples are listed in the table below.

Manual (AA1001) + CFX96			Magcore (AA1186) + CFX96			QiaSymphony(AA1439/192)+CFX96		
IU/ml	REPLICATES	POSITIVES	IU/ml	REPLICATES	POSITIVES	IU/ml	REPLICATES	POSITIVES
160	15	15	100	10	10	160	16	16
80	15	15	70	10	10	80	16	16
40	15	15	50	10	10	40	16	16
20	15	15	30	10	9	20	16	13
10	15	14	15	10	7	10	16	11
<b>LOD: 10.42 IU/ml</b>			<b>LOD: 32.86 IU/ml</b>			<b>LOD: 29.53 IU/ml</b>		

*Probit analysis at 95% done with Statplus 2009 .*

#### 2. Plasma

The analytical sensitivity was determined analyzing five scalar dilutions of a positive CMV plasma sample, calibrated with the First International CMV WHO Standard (NIBSC code 09/162). For serial dilutions a negative CMV human plasma was used. The analysis was carried out according to user manual, extraction included. The detection limits for CMV REAL TIME with plasma samples are listed in the table below.

Manual (AA1001) + CFX96			Magcore (AA1186) + CFX96			QiaSymphony (AA1440/96)+ CFX96		
IU/ml	REPLICATES	POSITIVES	IU/ml	REPLICATES	POSITIVES	IU/ml	REPLICATES	POSITIVES
100	17	17	1286	16	16	160	15	15
70	15	15	411	16	16	80	15	15
50	15	14	182	16	15	40	15	15
30	15	12	142	16	12	20	15	12
15	15	11	36,7	16	9	10	15	11
<b>LOD: 62.61 IU/ml</b>			<b>LOD: 43.44 IU/ml</b>			<b>LOD: 30.51 IU/ml</b>		

*Probit analysis at 95% done with Statplus 2009 .*

### Linear Range

To evaluate the linear range of the device, 8 dilution levels were prepared using real samples for dilutions from D to H; otherwise, for the points from A to C, the synthetic CMV DNA template (calibrated with the 1<sup>st</sup> WHO International Standard, NIBSC code 09/162), was diluted in CMV negative human blood and plasma. The CMV synthetic template was used due to the lack of clinical samples with high viral load.


The test was performed according to the instructions for use, extraction automatic included (AA1186, AA1439/192 and AA1440/96), using two different lots of reagents. Each dilution level has been tested in 10-15 replicates.


Dilution level	Concentration (IU/ml)	Concentration (log IU/ml)
A	5.00E+08	8.7
B	5.00E+07	7.7
C	5.00E+06	6.7
D	5.00E+05	5.7
E	5.00E+04	4.7
F	5.00E+03	3.7
G	5.00E+02	2.7
H	5.00E+01	1.7

Results:

**Blood**

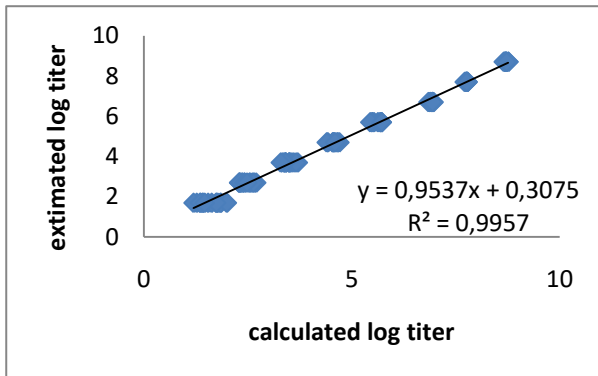**Plasma**

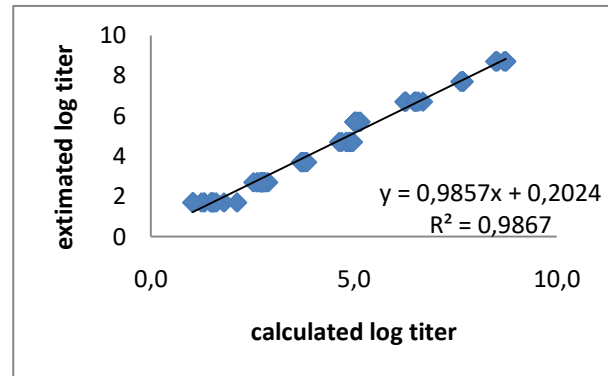Conclusions:

As shown in figure above, the device CMV REAL TIME showed a linear response from 5.00E+01 IU/ml to at least 5.00E+08 IU/ml, using acceptance criterion of  $\pm 0.3 \log_{10}$ .

**Diagnostic sensitivity**

The diagnostic sensitivity of the CMV REAL TIME kit was evaluated by analyzing 224 samples with a known positivity, extracted by NLM code AA1001, AA1186, AA1439/192 and AA1440/96 and amplified with the several Real Time instruments mentioned in the instruction for use.

In front of the results obtained, the diagnostic sensitivity of the device is 99.5%.

**Potentially cross reactive markers**

To evaluate the potential cross reactivity of CMV REAL TIME with other pathogens the in-silico analysis of primers and probes as well as the experimental one on real samples, coming from hospital routine were performed.

The in-silico study did not show cross reactivity with other pathogens, while the analysis of 34 samples positive for Herpesvirus, different from CMV, gave the following results:

Other pathogens	N° of samples	Positive with CMV Real-Time AA1489/48
EBV	6	0/6
BKV	6	0/6
HHV-6	5	0/5
HSV-1	4	0/4
HHV-8	1	0/1
HSV-2	4	0/4
VZV	8	0/8

⇒ All the samples gave the expected result. CMV REAL TIME did not show cross-reactivity with other Herpesvirus.

## ⇒ Precision

### Intra-assay reproducibility

Intra-assay precision was assessed using a total of 5 samples with different viral load. Two different operators analyzed each sample in triplicate using the same batch of reagents.

The results are shown in the table below:

<b>Viral load</b>	Sample1 4x10 <sup>4</sup> IU/ml	Sample2 1x10 <sup>4</sup> IU/ml	Sample3 4x10 <sup>4</sup> IU/ml	Sample4 1x10 <sup>3</sup> IU/ml	Sample5 9x10 <sup>2</sup> IU/ml
<b>CV%</b>	5.99	7.54	7.76	2.96	6.31

### Inter-assay reproducibility

Inter-assay precision was assessed using the same samples used for the intra-assay reproducibility. Two different operators tested 3 replicates in different runs and different days, using two different batches. The tests were performed according to the Instructions for use, starting from DNA extraction; in order to evaluate the precision for all the steps of the process.

The results are shown in the tables below:

<b>Viral load</b>	Sample1 4x10 <sup>4</sup> IU/ml	Sample2 1x10 <sup>4</sup> IU/ml	Sample3 4x10 <sup>4</sup> IU/ml	Sample4 1x10 <sup>3</sup> IU/ml	Sample5 9x10 <sup>2</sup> IU/ml
<b>CV%</b>	4.61	10.28	6.00	2.45	5.57

### **Whole system failure rate leading to false-negative results**

To evaluate the possibility of the presence of false negatives, 167 samples were analyzed, with viral load about 3 times higher than the detection limit for CMV.

The results are listed below and show a whole system failure rate <1%.

<b>Real Time Instrument used</b>	<b>CFX96</b>	<b>RGQ</b>	<b>AB 7500fast</b>
<b>Positive results</b>	48	54	64
<b>Negative Results for Positive Sample</b>	1	0	0

### **Cross-contamination**

⇒ The absence of cross-contamination between samples for the entire workflow was demonstrated by analysing 81 samples from different matrix in 3 different runs, alternating negative and viral loads samples.

All negative samples tested gave the expected negative result.

The CMV REAL TIME test did not show any cross-contamination between samples.


## INTERFERENT SUBSTANCES

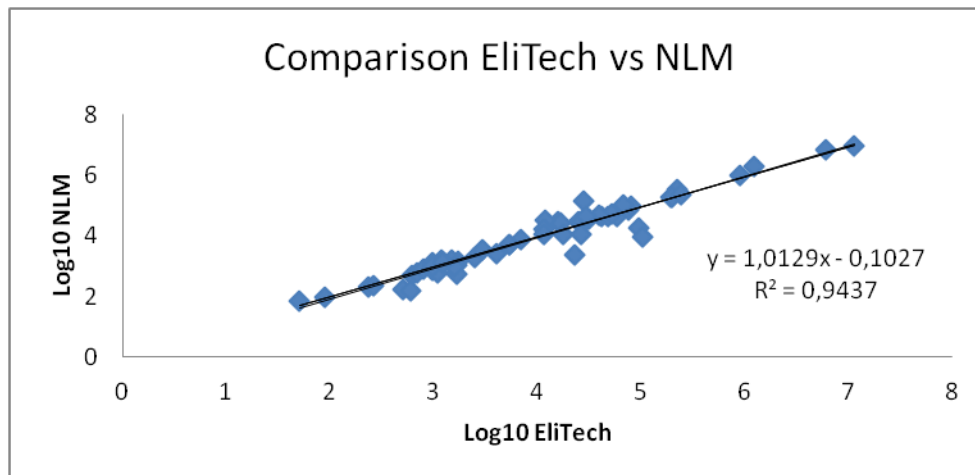
Data not available.

### Comparison with the CMV ELITE MGB® Kit in routine samples

The performance of the CMV REAL TIME device was compared to the **CMV ELITE MGB® Kit** system through the analysis of 63 samples reflecting the hospital routine conditions. The analysis was carried out in accordance with the user manual, including extraction.

The results revealed a good correlation between the two devices, as shown in the figure below:

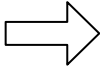						
Lot Number	Catalogue Number	Temperature Limitation	Use by	Single use	Positive Control	Consult Instructions for use
						
European Conformity	Manufacturer	<i>in vitro</i> diagnostic device	Sufficient for	No Template Control	Modified text compared with the previous version	