

BactoReal[®] Kit *Mycoplasma gallisepticum*

Manual

For use with the

- ABI PRISM[®] 7500 (Fast)
- Mx3005P[®]
- LightCycler[®] 480



For veterinary use only



DVEB00311, DVEB00313



100



DVEB00351, DVEB00353



50



ingenetix GmbH
 Arsenalstraße 11
 1030 Vienna, Austria
 T +43(0)1 36 1980 198
 F +43(0)1 36 1980 199
 office@ingenetix.com
 www.ingenetix.com

Index

1. Product description	3
2. Pathogen information	3
3. Principle of real-time PCR.....	3
4. General Precautions	3
5. Contents of the Kit.....	4
5.1. BactoReal® Kit <i>Mycoplasma gallisepticum</i> order no. DVEB00311 or DVEB00351.....	4
5.2. BactoReal® Kit <i>Mycoplasma gallisepticum</i> order no. DVEB00313 or DVEB00353.....	4
6. Additionally required materials and devices.....	4
7. Preparation of real-time PCR.....	5
7.1. Pipetting scheme	5
7.2. Programming of the temperature profile	5
8. Interpretation of PCR-data	6
9. Troubleshooting.....	8
10. Specifications	8
10.1. Analytical sensitivity.....	8
10.2. Analytical specificity.....	8
11. Annex – symbols.....	8

1. Product description

BactoReal® Kit *Mycoplasma gallisepticum* is a real-time PCR assay for detection of *M. gallisepticum* DNA. This test was developed and validated for the ABI PRISM® 7500 (Fast) instrument (Applied Biosystems), LightCycler® 480 (Roche) and Mx3005P® (Agilent), but is also suitable for other real-time PCR instruments. This test allows the rapid and sensitive detection of DNA of *M. gallisepticum* from samples purified from faecal samples, or swabs of tissues, exudates, diluted tissue homogenates, aspirates from the infraorbital sinuses or joint cavities, as well as from material from egg yolk or embryos. *Mycoplasma gallisepticum* DNA can be recovered efficiently from faecal samples using the QIAamp DNA Stool Mini Kit, and from swabs using the QIAamp DNA Mini Kit extraction methods, for example. *Mycoplasma imitans* is not detected by this test.

BactoReal® Kit *Mycoplasma gallisepticum* detects the *mgc2* gene (cytadhesin-like protein encoding gene) of *M. gallisepticum*. A probe-specific amplification-curve at 530 nm (FAM channel) indicates the amplification of *M. gallisepticum* specific DNA.

An internal positive control system for detection in VIC/HEX channel, (554 nm, order no. DVEB00311 or DVEB00351) or Cy5 channel (667 nm; order no. DVEB00313 or DVEB00353) excludes false-negative interpretation of results due to inhibition of real-time PCR (see 8. Interpretation of PCR-data).

When using PCR-platforms not validated by ingenetix, an evaluation of the multiplex-PCR is recommended. Please be aware that some PCR-platforms have to be calibrated with the corresponding dye before performing multiplex-PCR.

BactoReal®, MycoReal, ParoReal and ViroReal® Kits are optimized to run under the same thermal cycling conditions. RNA and DNA material can be analysed in one run.

2. Pathogen information

Avian mycoplasmosis is caused by several pathogenic mycoplasmas. *Mycoplasma gallisepticum* and *Mycoplasma synoviae* are the most important avian pathogens. In addition, *Mycoplasma meleagridis* and *Mycoplasma iowae* can also cause disease in poultry, but *Mycoplasma gallisepticum* and *Mycoplasma synoviae* are considered to be the main species of pathogenic mycoplasmas.

Mycoplasma gallisepticum and *Mycoplasma synoviae* cause chronic respiratory disease (CRD) and infectious synovitis of chickens and turkeys especially in the presence of stresses or other respiratory pathogens, causing economic losses due to decrease in growth rate, loss of weight and reduced egg production.

References:

Raviv, Z. and S.H. Kleven, 2009. The development of diagnostic Real-Time TaqMan PCRs for the four pathogenic avian mycoplasmas. Avian Dis., 53: 103-107.

3. Principle of real-time PCR

A specific DNA sequence of the pathogen genome is amplified and the generated PCR-product is detected by an oligonucleotide-probe labelled with a fluorescent dye. This technology allows for a sequence-specific detection of PCR amplicates.

4. General Precautions

The user should always pay attention to the following:

- Always include a negative control per PCR-run (water instead of sample).
- Optional: for valid interpretation of results, a negative control should be included during DNA-extraction (for example extraction of water instead of sample material), in order to exclude false-positive results due to contamination with *M. gallisepticum* DNA during extraction.
- Be careful when handling the positive control.
- Store and extract positive material (specimens, controls and amplicons) separately from all other reagents and add it to the reaction mix in a spatially separated workspace.
- Periodically decontaminate benches and devices.
- Use sterile pipette tips with filters.
- Thaw all components thoroughly at room temperature before starting an assay. When thawed, mix the components and centrifuge briefly.
- For MSDS, see www.ingenetix.com.

5. Contents of the Kit

5.1. BactoReal® Kit *Mycoplasma gallisepticum* order no. DVEB00311 or DVEB00351

Labelling	Content	Amount		Storage
		DVEB00311	DVEB00351	
<i>Mycoplasma gallisepticum</i> Assay Mix (green cap)	Primer and probe (FAM) for detection of <i>M. gallisepticum</i>	2 x 50 µl	1 x 50 µl	-20°C
CR-1 Assay Mix (yellow cap)	Primer, probe (VIC/HEX) and target for detection of IPC	2 x 50 µl	1 x 50 µl	-20°C
<i>Mycoplasma gallisepticum</i> Positive Control (red cap)	Control-DNA (approx. 10,000 target copies/µl)	1 x 25 µl	1 x 25 µl	-20°C
DNA Reaction Mix (white cap) [#]	Reaction Mix	2 x 500 µl	1 x 500 µl	-20°C until first use, then at +4°C
Water (blue cap)	Water	1 x 1000 µl	1 x 1000 µl	-20°C to +4°C

[#]DNA Reaction Mix contains uracil-N glycosylase (UNG)

5.2. BactoReal® Kit *Mycoplasma gallisepticum* order no. DVEB00313 or DVEB00353

Labelling	Content	Amount		Storage
		DVEB00313	DVEB00353	
<i>Mycoplasma gallisepticum</i> Assay Mix (green cap)	Primer and probe (FAM) for detection of <i>M. gallisepticum</i>	2 x 50 µl	1 x 50 µl	-20°C
CR-3 Assay Mix (yellow cap)	Primer, probe (Cy5) and target for detection of IPC	2 x 50 µl	1 x 50 µl	-20°C
<i>Mycoplasma gallisepticum</i> Positive Control (red cap)	Control-DNA (approx. 10,000 target copies/µl)	1 x 25 µl	1 x 25 µl	-20°C
DNA Reaction Mix (white cap) [#]	Reaction Mix	2 x 500 µl	1 x 500 µl	-20°C until first use, then at +4°C
Water (blue cap)	Water	1 x 1000 µl	1 x 1000 µl	-20°C to +4°C

[#]DNA Reaction Mix contains uracil-N glycosylase (UNG)

The components of BactoReal® Kit *Mycoplasma gallisepticum* are stable until the expiry date stated on the label. Repeated thawing and freezing should be avoided. Please protect kit components from light.

6. Additionally required materials and devices

- Reagents and devices for DNA-extraction
- PCR-grade water
- Disposable powder-free gloves
- Pipettes (adjustable)
- Sterile pipette tips with filters
- Vortex mixer
- Desktop centrifuge with rotor for 2 ml reaction tubes
- Real-time PCR instrument which is able to detect and differentiate fluorescence in FAM and VIC/HEX or Cy5 channel
- Appropriate 96 well reaction plates or reaction tubes with corresponding (optical) closing material

7. Preparation of real-time PCR

Please make sure that at least one negative control (water, blue cap), as well as one positive control (red cap) and one extraction negative control (optional, recommended) are included per PCR run.

Ingenetix highly recommends performing PCR analyses in duplicates, which increases the probability of detection of the pathogen and facilitates interpretation of results.

7.1. Pipetting scheme

		Per sample
Preparation of Master Mix (mix well)	Water*	3.0 µl
	DNA Reaction Mix (2x)	10.0 µl
	<i>Mycoplasma gallisepticum</i> Assay Mix	1.0 µl
	CR Assay Mix	1.0 µl
	Total volume Master Mix	15.0 µl
Preparation of PCR	Master Mix	15.0 µl
	Sample*	5.0 µl
	Total volume	20.0 µl

*1-8 µl of the sample can be used. When using an amount other than 5 µl of the sample, the amount of H₂O has to be changed accordingly.

Positive Control: As positive control use 1 µl of the *Mycoplasma gallisepticum* Positive Control + 4 µl H₂O.

Optional: a 1:10 dilution of the positive control can be used and defined as second standard value (approx. 1000 target copies/µl).

7.2. Programming of the temperature profile

Please find further information on programming the real-time PCR instrument in the respective operator's manual. Please be aware that some PCR-platforms have to be calibrated with the corresponding dye before performing multiplex-PCR.

Select dyes: FAM-TAMRA for detection of *Mycoplasma gallisepticum*

Cy5-NONE (CR-3 Assay Mix) or VIC-TAMRA (CR-1 Assay Mix) for detection of IPC

Select reference dye (passive reference): ROX

Sample Volume: 20 µl

Temperature Profile:

Program 1	Program 2	Program 3
Cycles: 1 Analysis: None	Cycles: 1 Analysis: None	Cycles: 45 Analysis: Quantification Acquisition at 60°
50°C 2 min*	95°C 20 sec	95°C 5 sec
		60°C 1 min

For ABI PRISM® 7500:

Ramp speed: Without "fast cycling" parameter

For LightCycler® 480 instrument:

Detection format: 2 Color Hydrolysis Probe
(dyes see above)

***Note:** If viral RNA should be also detected in the same PCR run, program 1 has to be prolonged to 15 min at 50°C. This temperature profile can be used for all BactoReal®, MycoReal®, ParoReal® and ViroReal® kits for the detection of DNA or RNA.

8. Interpretation of PCR-data

Examples for interpretation of positive reactions are shown in the amplification plots below.

For a valid interpretation, the following criteria must be fulfilled:

	Ct/Cp (FAM channel) <i>M. gallisepticum</i> target	Ct/Cp IPC target	Interpretation
Negative control	Negative	36.0 ± 2	Valid
Positive control (undiluted, 1 µl/PCR)	28.0-31.0	36.0 ± 2	Valid
Or: positive control (1:10 diluted, 1 µl/PCR)	31.0-34.0	36.0 ± 2	Valid
Extraction negative control (optional)	Negative	36.0 ± 2	Valid
Negative sample	Negative	36.0 ± 2	Valid
Positive sample	Positive	Positive/Negative	Valid

For analysis of PCR data please proceed as follows:

For analysis of PCR results gained with BactoReal® Kit *Mycoplasma gallisepticum* please select fluorescence display options FAM channel for the *M. gallisepticum* target and VIC/HEX channel (order no. DVEB00311, DVEB00351) or Cy5 channel (order no. DVEB00313, DVEB00353) for the internal positive control target. Samples with a positive Cp or Ct-value are considered positive. Please also check the presence of amplification-curves manually.

Once the analysis is completed, the following results are possible:

1. Signal in FAM channel:

→ DNA of *M. gallisepticum* was amplified. The sample has to be interpreted as positive.

M. gallisepticum DNA can lead to a reduced or absent fluorescence signal of the internal positive control (competition of PCR).

2. No signal in FAM channel:

→ No *M. gallisepticum* DNA is detectable in the sample. The sample has to be interpreted as negative. An inhibition of PCR cannot be excluded.


2a. No signal in FAM channel but signal of the internal positive control:

→ No *M. gallisepticum* DNA is detectable in the sample. The sample has to be interpreted as negative. The positive signal of the internal positive control assay excludes a putative PCR inhibition.

2b. No signals in FAM channel and no signal with internal positive control:

→ No interpretation statement can be made.

Information about possible sources of error and their solution can be found in 9. Troubleshooting.

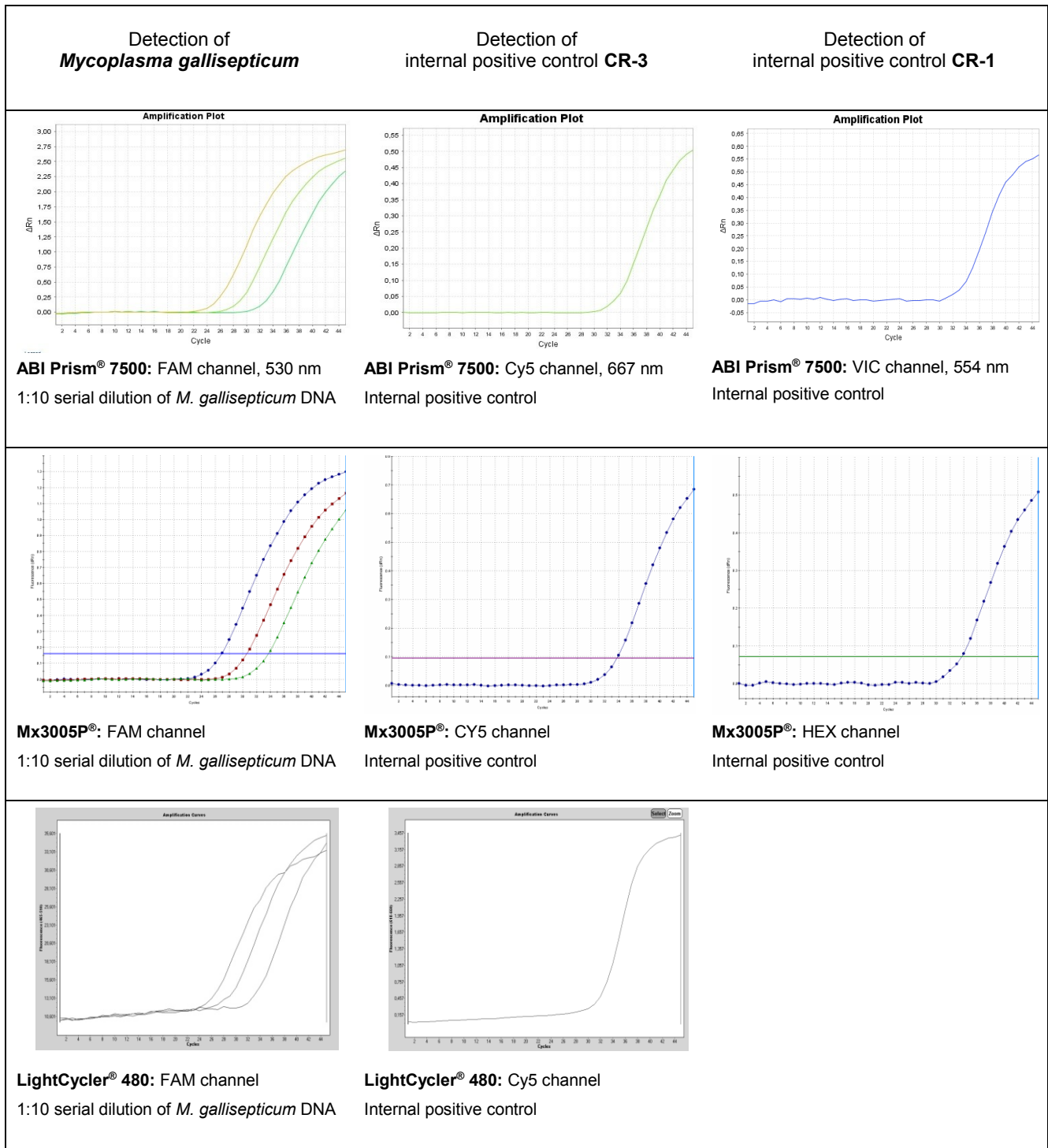

9. Troubleshooting

1. No *M. gallisepticum* specific signal with positive control:

- Incorrect programming of the temperature profile of the real-time PCR instrument.
→ Compare the temperature profile with the protocol (see 7. Preparation of real-time PCR).
- Incorrect configuration of the PCR reaction.
→ Check your work steps (see 7. Preparation of real-time PCR) and repeat the PCR, if necessary.

2. No signal with the internal positive control and no *M. gallisepticum* specific signal with the sample:

- The PCR reaction was inhibited. No interpretation can be made.
→ Make sure that you use a recommended method for DNA isolation and stick closely to the manufacturer's instructions.
→ If no operating mistakes during extractions can be retraced, it is recommended to repeat the PCR with lower amounts of DNA-eluate (1/5 or 1/10 of sample volume + the adequate amount of H₂O).
- Incorrect PCR conditions.
→ Check the PCR conditions and repeat the PCR, if necessary.

3. *M. gallisepticum* specific signal with the negative control:

- A contamination occurred during preparation of the PCR.
→ Repeat the PCR with new reagents in replicates.
→ Strictly pipette the positive controls at last.
→ Make sure that work space and instruments are decontaminated at regular intervals.

4. *M. gallisepticum* specific signal with the negative control of DNA-extraction:

- A contamination occurred during extraction.
→ Repeat the extraction and PCR using new reagents.
→ Make sure that work space and instruments are decontaminated at regular intervals.

10. Specifications

BactoReal® Kit *Mycoplasma gallisepticum* was evaluated with the ABI PRISM® 7500 (Fast) instrument (Life Technologies), with the LightCycler® 480 (Roche) and the Mx3005P® (Agilent). For further validation data please contact ingenetix.

10.1. Analytical sensitivity

The analytical sensitivity is 10 target copies/PCR reaction. The limit of detection (LoD₉₅ = smallest number of copies of target DNA which can be detected in 95% of cases) is 27 target copies/reaction and was determined by several replicates around the detection limit.

10.2. Analytical specificity

The specificity is ensured by the selection of highly specific primers and probes. The primers and probes were checked for possible homologies to currently published sequences by sequence comparison analyses. This also validated the detection of so far known *M. gallisepticum* strains.

Specificity was tested on isolates of *M. gallisepticum*, *M. meleagridis* and *M. synoviae*. No cross reactions were observed. *Mycoplasma imitans*, a *Mycoplasma* species that is serologically related to *M. gallisepticum* and that presents the same biochemical properties is not detected by this test. Eight field samples were tested and correctly analysed.

11. Annex – symbols



Batch code



Catalogue number



Contains sufficient for <n> tests



Use by



Manufactured by



Store at