


biotical

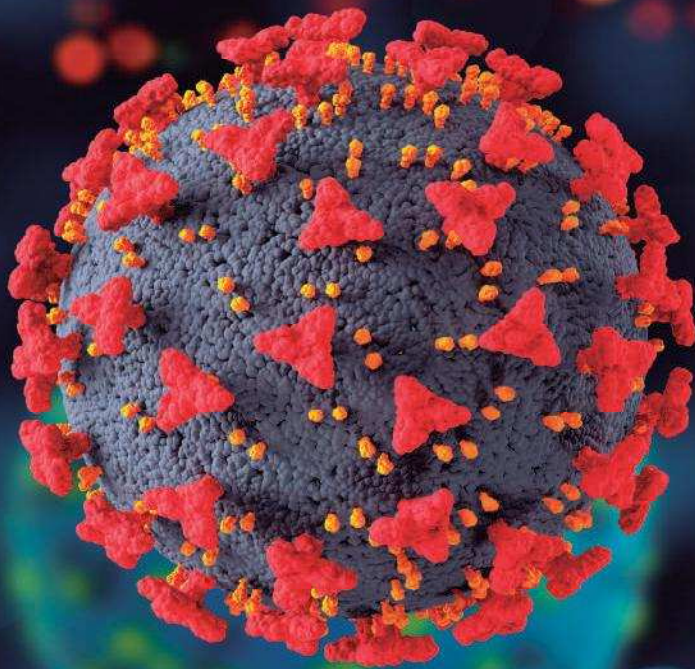

IN VITRO DIAGNOSTIC DEVICES
ready to test

ANTIGEN RAPID TEST

biocore diagnostics GmbH
Louisen Str. 137
D-61348 Bad Homburg,
Germany


Tel: +49-6172-17102-16
Fax: +49-6172-17102-29
E-mail: info@biocore-diagnostics.de
www.biocore-diagnostics.de



Commerzbank Bad Homburg
Bank Code: 500 400 00
WIFT No.: COBADEFFxxx,
Acc. No.: 3467560
IBAN: DE13 5004 0000 0346 756000

Tax No. 00322950108
VAT. ID No. DE 160003733
Trade Register: HRB 11589
Managing Directors:
Dr. Gert Sokolowski
Mr. Adnan Dawoud

biotical health

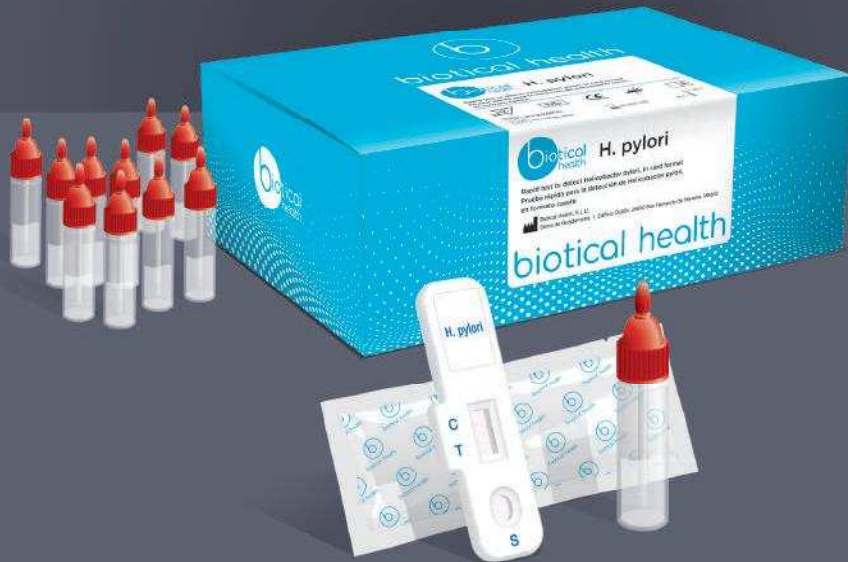

devices ready to test

- All our tests are manufactured and developed in Spain (including monoclonal antibodies)
- No equipment required
- We can supply all our references with positive and negative controls, for extra quality assurance
- Latex particle technology to obtain a different colour in the control and test line, to facilitate the interpretation of the result
- No cold chain required
- The same sample can be used for several test in the same category

TEST CATALOGUE


	TEST REFERENCE	POSITIVE CONTROL REF.	
RESPIRATORY PATHOGENS	NASAL SAMPLE		
	RSV	RTB25RVC	CPRT30RV
	ADENO RESPIRATORY	RTB25APC	CPRT30AP
	RSV+ADENO RESPIRATORY	RTB25VRC	-
	INFLUENZA A + B	RTB25ABC	CPRT30AB
	SARS-COV-2 AG	RTB25COV2C	CPRT30COV2
	OROPHARYNGEAL SAMPLE		
	STREP A	RTB25SC	CPRT30S
	URINE SAMPLE		
	S. PNEUMONIAE	RTB25SNCN	CPRT30SN
LEGIONELLA	RTB25LPCN	CPRT30LP	
GASTROINTESTINAL PATHOGENS	VIRUSES		
	ROTAVIRUS	RTB25R	CPRT30R
	ADENO INTESTINAL	RTB25AD	CPRT30AD
	ROTAVIRUS + ADENOVIRUS	RTB25RA	-
	ASTROVIRUS	RTB25AT	CPRT30AT
	NOROVIRUS GEN I + II	RTB25NV	CPRT30NG NOR. GEN I
			CPRT30NP NOR. GEN II
	ENTEROVIRUS	RTB25EV	CPRT30EV
	BACTERIA		
	E.COLI	RTB25E	CPRT30E
	CAMPYLOBACTER	RTB25CA	CPRT30CA
	SALMONELLA	RTB25SA	CPRT30SA
	SALMONELLA TYPHI	RTB25ST	CPRT30ST
	SALMONELLA PARATYPHI	RTB25SP	CPRT30SP
	H. PYLORI	RTB25PG	CPRT30PG
	C. DIFFICILE GDH	RTB25GD	CPRT30GD
	C. DIFFICILE TOXINS A + B	RTB25CD	CPRT30TA C. D. TOXIN A
			CPRT30TB C. D. TOXIN B
	YERSINIA ENTEROCOLITICA 0:3	RTB25Y3	CPRT30Y3
	YERSINIA ENTEROCOLITICA 0:9	RTB25Y9	CPRT30Y9
	PARASITE		
	GIARDIA	RTB25G	CPRT30G
	CRYPTO	RTB25K	CPRT30K
	CRYPTO + GIARDIA	RTB25KG	-
	ENTAMOEBIA	RTB25EH	CPRT30EH
TUMORS AND INFLAMMATORY MARKERS	FOB	RTB25F	CPRT30F
	TRANSFERRIN	RTB25TF	CPRT30TF
	FOB + TRANSFERRIN	RTB25FT	-
	CALPROTECTIN	RTB25CP	CPRT30CP
	LACTOFERRIN	RTB25LT	CPRT30LT
	CALPROTECTIN + LACTOFERRIN	RTB25CL	-
ALL TEST NEGATIVE CONTROL	CNRT30		


Respiratory pathogens


biocore diagnostics

Our Experience Is Your Advantage

Influenza A+B test procedure


- 

1 Add 15 drops of the Reagent B in a testing tube.
- 

2 Put the swab into the tube and mix the solution by rotating the swab forcefully against the side of the tube at least 1 minute.
- 

3 Dispense exactly 4 drops from the testing tube into the "S" window.

Read the result at 10 minutes

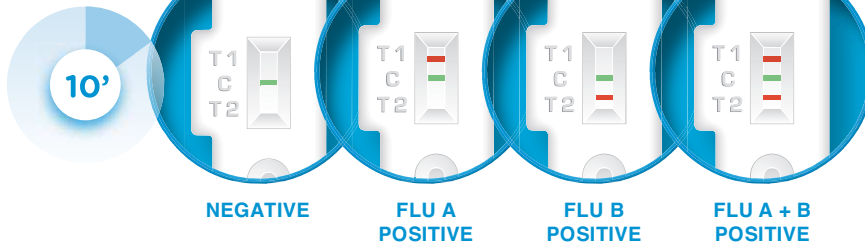


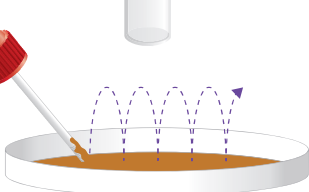
Gastrointestinal pathogens


E. coli test procedure

Unscrew the cap and use the stick to pick up the sample.

Cut the end of the cap and dispense exactly 4 drops into the "S" window.


- 

1
- 

2 Close the tube with the diluent and stool sample. Shake the tube in order to assure good sample dispersion.
- 

3

Read the result at 10 minutes

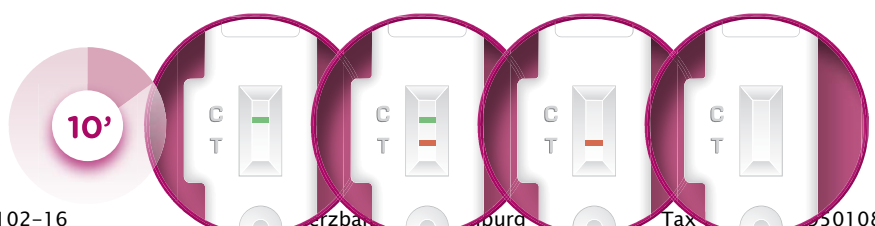


biocore diagnostics GmbH
Louisen Str. 137
D-61348 Bad Homburg,
Germany

Tel: +49-6172-17102-16
Fax: +49-6172-17102-29
E-mail: info@biocore-diagnostics.de
www.biocore-diagnostics.de

Bank Code: 500 400 00
SWIFT No.: COBADEFFxxx,
Acc. No.: 3467560

Tax No.: 350108
VAT. ID No. DE 160003733
Trade Register: HRB 11589
Managing Directors:


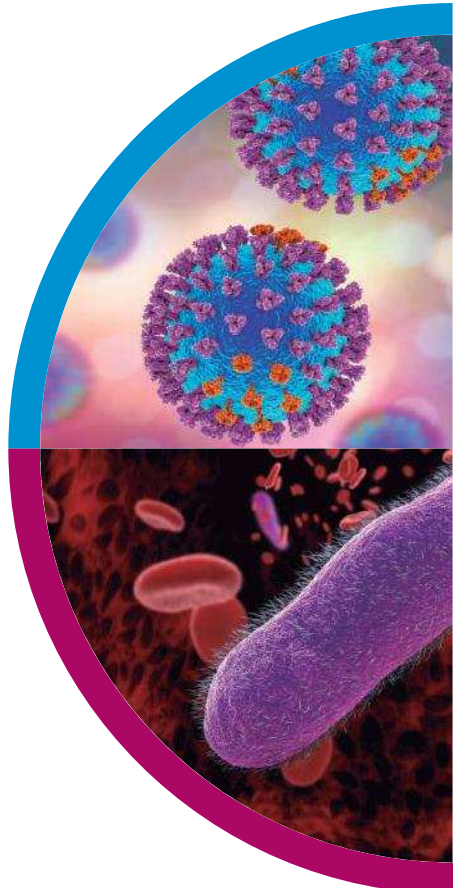


Tumor and inflammatory markers biocore diagnostics

Our Experience Is Your Advantage


Calprotectin test procedure

Unscrew the cap and use the stick to pick up the sample.

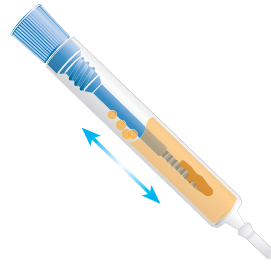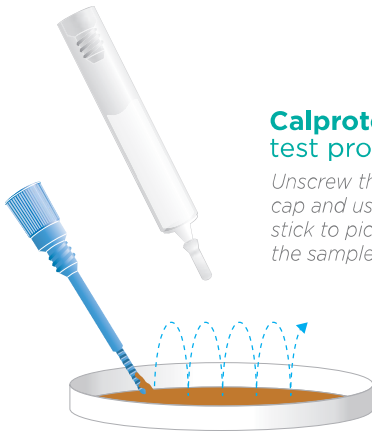

1 2 3

Close the tube with the diluent and stool sample. Shake the tube in order to assure good sample dispersion.

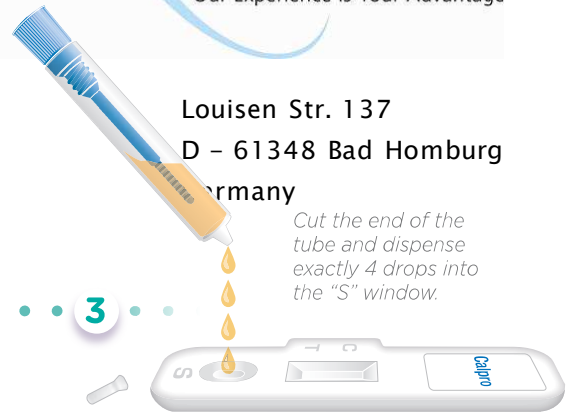

Louisen Str. 137
D - 61348 Bad Homburg
Germany

Cut the end of the tube and dispense exactly 4 drops into the "S" window.

10'

Read the result at 10 minutes

NEGATIVE

POSITIVE

INVALID

INVALID

Controls (+&-)

+

Results for the Strep. A positive control

POSITIVE CONTROL

INVALID

INVALID

INVALID

-

Results for the negative control

NEGATIVE CONTROL

INVALID

INVALID

INVALID



biocore diagnostics GmbH
Louisen Str. 137
D-61348 Bad Homburg,
Germany

Tel: +49-6172-17102-16
Fax: +49-6172-17102-29
E-mail: info@biocore-diagnostics.de
www.biocore-diagnostics.de

SWIFT NO.: LUBADEFXXX,
Acc. No.: 3467560

Trade Register: HKB 11
Managing Directors: